

(19) World Intellectual Property Organization  
International Bureau




(43) International Publication Date  
11 December 2003 (11.12.2003)

PCT

(10) International Publication Number  
**WO 03/103010 A1**

- (51) International Patent Classification<sup>7</sup>: **H01J 49/42**, 49/40 (74) Agent: **HEISLER, Bradley, P.**; Heisler & Associates, Suite 300, 3017 Douglas Blvd., Roseville, CA 95661 (US).
- (21) International Application Number: PCT/US03/15718 (81) Designated States (*national*): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NI, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW.
- (22) International Filing Date: 15 May 2003 (15.05.2003)
- (25) Filing Language: English
- (26) Publication Language: English
- (30) Priority Data: 10/159,222 31 May 2002 (31.05.2002) US (84) Designated States (*regional*): ARIPO patent (GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HU, IE, IT, LU, MC, NL, PT, RO, SE, SI, SK, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).
- (71) Applicant (*for all designated States except US*): **MICHROM BIORESOURCES, INC.** [US/US]; 1945 Industrial Drive, Auburn, CA 95603 (US).
- (71) Applicants and (72) Inventors: **NUGENT, Kerry, D.** [US/US]; 10560 Canyon Creek Lane, Penn Valley, CA 95946 (US). **WANG, Houle** [US/US]; 12461 Morning Mist Lane, Auburn, CA 95602 (US).
- Published:**  
— with international search report
- For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.*

(54) Title: TWO-DIMENSIONAL TANDEM MASS SPECTROMETRY

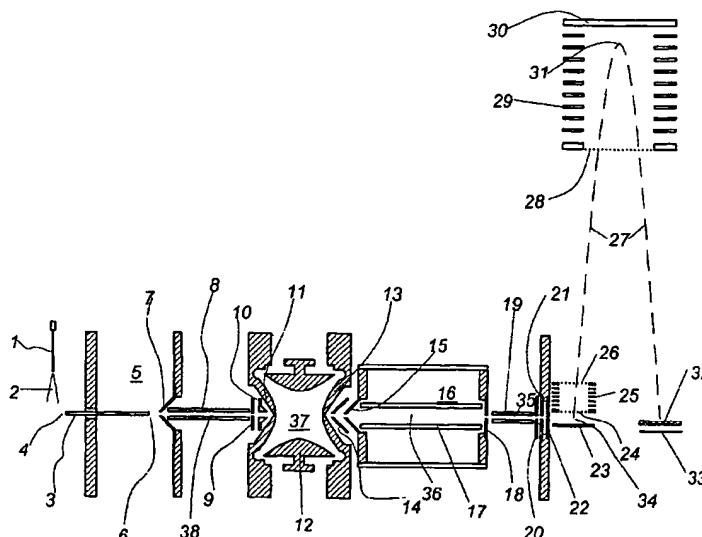

(57) Abstract: A tandem mass spectrometer is provided including two mass analyzers with an ion fragmentation device interposed between the two mass analyzers. The first mass analyzer is a non-destructive mass analyzer, such as an ion trap, to initially collect and hold parent ions and sequentially release parent ions of known mass to charge ratio. The released parent ions pass through the fragmentation device, such as a collision cell, where the parent ions are fragmented into daughter ions. These daughter ions then pass on to the second mass analyzer. The second mass analyzer is of a high speed full spectrum type, such as a time of flight analyzer, so that a full spectrum of mass data is provided for the daughter ions, to go with parent ion mass spectrum data from the first mass analyzer.

WO 03/103010 A1

**TWO-DIMENSIONAL TANDEM MASS SPECTROMETRY**

5

**Technical Field**

10 The present invention relates to mass spectrometry apparatuses and methods for obtaining data which identify the mass to charge ratio of various parent ions in a sample as well as mass to charge ratio of daughter ions produced by fragmentation of the parent ions in the sample, such as to determine structural information about the parent ions, and to derive other information about relationships between the parent ions and daughter ions. More particularly, this invention relates to mass spectrometry systems which include tandem mass analyzers separated by an ion  
15 fragmentation cell to obtain multi-dimensional data about the parent ions and daughter ions of the sample.

**Background Art**

20

In simple mass spectrometers (MS), sample ions are formed in an ion source, such as by Electron Impact (EI), or by Atmosphere Pressure Ionization (API). The ions then pass through a mass analyzer, such as a quadrupole or time of flight device (TOF), for detection. The detected ions can be molecular ions (parent ions), fragment ions (daughter ions) of the molecular ions, or  
25 fragment ions of other daughter ions.

Quadrupole mass analyzers and magnetic sector mass analyzers, are mass filter type mass analyzers that allow only ions with specific mass/charge ratios ( $m/z$ ) to pass through. Other ions are discarded during the scan. These type of mass analyzer is not non-destructive. This type of mass analyzer is thus not particularly effective for a full mass scan (also called full spectrum scan)  
30 where multiple ions of different  $m/z$  in a sample are to be detected and/or measured. Ion trap mass analyzers can trap ions and then analyze them sequentially based on the Fourier Transform Ion Cyclotron Resonance (FT-ICR)  $m/z$ . Mass analyzers can obtain similar full spectrum data, but in a different fashion by first measuring all of the ions and then performing a fourier transform analysis to measure the different ions in the sample. Therefore, the duty cycle and effectiveness of these  
35 types of non-destructive mass analyzers for full mass scans is higher than for mass filter type instruments. Time of flight mass analyzers sort ions based on flight time from an accelerator region to a detector spaced from the accelerator region. TOF mass analyzers can detect all ions, no matter what their mass to charge ratios are, and so they have very good sensitivity for a full mass scan spectrum.

Ion fragmentation mass spectrometers have been developed, characterized by having two or multiple sequential stages of mass analysis and an intermediate fragmentation region where parent ions from the first stage are fragmented into daughter ions for the second stage. Hence, these are generally termed "tandem" or "MS/MS" instruments. In such tandem mass spectrometers, sample ions are produced in an ion source, and the first stage of mass analysis analyzes selected parent ions of particular mass or  $m/z$  with a mass filter type mass analyzer. Then, some of the selected parent ions are fragmented or otherwise caused to dissociate, such as by metastable decomposition, collision induced dissociation (CID), or collisionally activated dissociation (CAD), to produce the daughter ions. Finally, the second stage of mass analysis sorts the daughter ions according to mass or  $m/z$ .

There are two styles of instruments in terms of "tandem" mass spectrometers, "tandem in space" and "tandem in time." Tandem in space mass spectrometers, such as triple quadrupoles and quadrupole-time of flight (Q-TOF) devices, have two mass analyzers, one for parent ion selection and one for daughter ion detection and/or measurement. Two mass analyzers are separated by a fragmentation device. Tandem in time instruments, on the other hand, have one mass analyzer that analyses both parent ions and daughter ions, but sequentially in time. Ion trap and FT-ICR are two most common mass spectrometers that have tandem in time MS/MS. The parent ions first are selected in the analyzer cell then fragmented. Often fragmentation takes place inside the analyzer. Then the daughter ions are analyzed in the same cell. Alternatively it is known to analyze the daughter ions in a downstream analyzer, such as a TOF analyzer.

Several MS/MS scan types are used based on the relationship between the parent ions and the daughter ions. "Daughter scan" is a method that involves a full scan of daughter ions while the parent ion from which the daughter ions originate is pre-selected and fixed. This method is useful if an analyst knows the molecular weight of the parent ion and wants to know structural information about the parent ion. For instance, two distinct parent ions of similar molecular weight, but different structure can be differentiated by what daughter ions they typically fragment into. The data dependent daughter scan is often used when combined with liquid chromatographs (LC-MS/MS). The mass spectrometer automatically selects a parent ion peak based on previous scans and the peak intensity, charge state and other considerations. The mass analyzer then makes a full scan of the daughter ions resulting from fragmentation of the parent ion of interest.

"Parent ion scan," also known as "precursor scan," is a method that has a fixed daughter ion selection for the second analysis stage, while using the first stage to scan all of the pre-fragmentation parent ions in the sample. Only those molecules/compounds in the sample are detected which produce a specific daughter ion when fragmented. If both parent ion selection and daughter ion selection are fixed, an analyst will get selected reaction monitoring (SRM). SRM has the best selectivity and good signal to noise ratio for quantitation.

"Neutral loss scan" is a method that shows all parent ions that lose a particular mass during fragmentation. The second stage mass analyzer scans the ions together with the first stage mass

analyzer but with a certain offset. Neutral loss scans are used for screening experiments where a group of compounds all give the same loss.

5 Magnetic and electrostatic sector (together referred to as "sector") mass analyzers have relatively slow scan speed, so sector based MS/MS instruments including sector-sector, sector-quadrupole and sector-TOF are normally good for daughter scans which don't need high speed scanning of parent ions in the first stage. Tandem in time instruments select the parent ion first, then fragment and scan the daughter ions later. Normally this type of instrument can only perform full mass scan of the daughter ions.

10 Time of flight mass analyzers are known to have a number of advantages, including fast scanning rate, higher sensitivity, relatively high resolution and good mass accuracy. Q-TOF is a MS/MS instrument that combines quadrupole and TOF analyzers. It gives very good mass accuracy and sensitivity on full mass daughter scans but only filters a chosen parent ion with other parent ions being lost.

15 Triple quadrupole mass spectrometers can do all of the above scans. However, since both the first and second stages of mass analysis are of the mass filter type, triple quadrupole systems are generally less effective than ion trap for full scan MS/MS, and less accurate and sensitive than Q-TOF.

20 To solve modern analytical problems an analyst often needs to use more than one MS/MS scan method. For LC-MS/MS the parent ions duration time is limited because additional peaks elute from the LC device in a specified time period. Normally there is not enough time to do different types of scans in a single LC run. It is also not unusual that several parent ions co-elute at the same time. In many cases, data dependent scans do not have enough time to fully analyze all parent ions.

25 A combined sector and TOF mass spectrometer is described in Enke et al U.S. Pat. No. 4,472,631. In Enke's method, a collision cell is placed before a magnetic sector. A pulsed ion source is also used, so that the flight time of the ion can be measured. The time resolution is used for parent ion information while a spatial resolution from a sector is used to give daughter ion information. By using a digital computer, a partial two dimensional spectrum of the selected parent ion and daughter ions can be reconstructed.

30 In Enke's invention, two spatial scan methods are described. One uses a fixed slit before the ion detector. Different daughter ion spectrums can be obtained by scanning magnetic field strength on the sector. For this method, only daughter ions with a particular  $m/z$  can be detected at a time. Daughter ions with a  $m/z$  other than this particular range of  $m/z$  will be thrown away. Less than 1% of all possible useful information can be obtained by the Enke device. This device is thus not effective to obtain highly sensitive full scan daughter ion spectrums.

A design using a multi-channel spatial array detector is also described by Enke. With this design, magnetic field strength within the magnetic sector is not scanned during operation. Rather, a micro-channel array, positioned at the focal plane of the magnetic sector, simultaneously detects

and individually resolves ion currents from a plurality of ion paths by use of individual micro-channels. The individual outputs of the micro-channel array are connected through amplifiers to individual time array detectors, connected to a digital computer. This method provides much better detection efficiency with a high duty cycle, but the spatial resolution is limited by the number of  
5 detector arrays and the size of the instrument. For a high resolution measurement, thousands of detector elements and associated electronics would be needed.

### Disclosure of Invention

10

Parent ions are first separated by a relatively slow, non-destructive scan device, for example, an ion trap. These parent ions are collected within the ion trap and then selectively released into a fragmentation device, such as a collision cell external to the first analyzer. Parent ion information is determined based on the time that individual parent ions are released from the ion trap or other first  
15 mass analyzer. The fragmentation devices sequentially fragment the parent ions into daughter ions. Then each daughter ion is analyzed by a fast scan analyzer, for example, a time of flight (TOF) mass analyzer.

In TOF scan, all ions from the same scan are originally from parent ions having the same mass/charge ratio ( $m/z$ ). In a certain range, all ions will be fragmented and scanned by TOF scans.  
20 A complete two-dimensional MS/MS map can be obtained after a single ion trap scan. A full scan MS spectrum can also be reconstructed by plotting total ion counts for each TOF scan.

Different MS/MS scans such as daughter scan, parent scan, neutral loss scan and selected reaction monitoring are all subsets of this complete 2-D MS/MS map.

During the MS/MS scan, unlike ion filter type instruments, no unnecessary ion loss occurs. A  
25 multi-pole ion guide with an electric ion gate prior to the ion trap can also act as an ion reservoir during the scan. Therefore, a theoretical 100% efficiency can be achieved.

### Brief Description of Drawings

30

FIG. 1 is a block diagram of a two-dimensional ion trap-TOF tandem mass spectrometer with an external collision cell.

FIG. 2 is a block diagram of two-dimensional ion trap-TOF tandem mass spectrometer with an external infrared multi-photon dissociation (IRMPD) cell.

35 FIG. 3 is a timing diagram that shows the correlation between the first stage analyzer and the second stage analyzer of the tandem mass spectrometer.

FIG. 4 is a three-dimensional graphical MS/MS map of a mixture of five different angiotensins shown simulating one example of how the MS/MS map of this invention would appear.

FIG. 5 is a two-dimensional plot of the MS/MS map of FIG. 4, viewed from above, showing the different subsets of MS/MS scans.

FIG. 6 is a daughter ion scan subset of the MS/MS map of FIG. 5 for a single parent ion ( $m/z=884$ ) simulating how such a daughter ion scan would appear using the subset two-dimensional MS/MS of this invention.

FIG. 7 is a parent ion scan subset of the MS/MS map of FIG. 5 for a single daughter ion ( $m/z=610$ ) simulating how such a parent ion scan would appear using the two-dimensional MS/MS of this invention.

FIG. 8 is a neutral loss scan subset of the MS/MS map of FIG. 5 simulating how such a scan would appear using the two-dimensional MS/MS of this invention.

FIG. 9 is a neutral loss two-dimensional map representing the X-axis in terms of amount of neutral loss.

FIG. 10 is a re-constructed full scan first stage MS spectrum of all of the parent ions, simulating how such a scan would appear using the two dimensional MS/MS of this invention.

### Best Modes for Carrying Out the Invention

Referring to the drawings, wherein like reference numerals represent like parts throughout the various drawing figures, FIG. 1 depicts a tandem mass spectrometer featuring an ion trap as a first mass analyzer and a time of flight device as a second mass analyzer according to a preferred embodiment of this invention. The two mass analyzers are separated by a fragmentation cell. In FIG. 2, a variation on the tandem mass spectrometer of FIG. 1 is shown where an infrared laser is included as part of the fragmenter between the two mass analyzers.

In essence, and with particular reference to FIG. 1, a sample is typically first ionized and then fed into an ion guide within a vacuum region leading the ions of the sample into the ion trap or other first stage mass analyzer. The ion trap thus contains one or more species of parent ions therein. As a voltage of the ion trap is increased, ions of different mass/charge ratio ( $m/z$ ) are sequentially released from the ion trap with such release detected so that a mass/charge ratio for the ions being released is determined. The parent ions released from the ion trap are then passed through a fragmentation cell, where various different fragmentation methodologies can be utilized to divide the parent ions passing therethrough into daughter ions. These daughter ions are then passed on to a second stage mass analyzer preferably in the form of a time of flight (TOF) mass analyzer. The TOF mass analyzer accelerates the daughter ions and then measures an amount of time from ion acceleration until impacting a detector. This time is correlated with the mass/charge ratio of the daughter ions.

Data collection, preferably in the form of a digital computer, is coupled to the ion trap mass analyzer and the TOF mass analyzer so that two dimensional data representative of the mass/charge

ratios of both the parent ions and the daughter ions (i.e. FIG. 5), as well as the relative abundance potentially forming a third dimension (FIG. 4), can be plotted a variety of different ways.

More specifically, and with particular reference to FIG. 1, details of the tandem mass spectrometer according to a preferred embodiment of this invention, is described. In FIG. 1 a liquid sample is ionized, such as via electrospray, by applying a high voltage between an electron spray ionization (ESI) needle 1 and the end of sample inlet capillary 4. Charged droplets and/or gaseous phase ions pass through the sample capillary 3 and enter into the low vacuum region 5 which is pumped by a roughing pump to about 1 mbar. Most of the air, moisture and neutral solvent molecules are pumped away in this stage. A cone shaped skimmer 7 allows ions through to the next stage. Preferably, a RF only multi-pole ion guide 8 is placed in the next pump region. The pressure in this region is between 0.01 to 0.001 mbar. In such pressure, ions will undergo collisional cooling 38. An electrostatic lens 9, 10 is preferably provided to further focus the ion beam. The above ion source details are typical, but any technique for delivering sample ions to the first stage analyzer of this invention can be similarly utilized.

The ion trap itself can be of either a three dimensional variety or configured as a linear ion trap. Preferably, the ion trap is of the three dimensional type and includes two end cap electrodes 11, 13 and a ring electrode 12 which together form an electric field to trap the parent ion therein. Ions pass through an ion trap inlet, typically in the form of a hole in the end cap electrode 11 and are first trapped in center region 37. These parent ions in the sample are then sequentially released through an ion trap outlet, typically in the form of an exit hole in end cap 13, based on their mass charge ratio  $m/z$ . Before the parent ions enter an entrance of the collision cell 16, the kinetic energy of the parent ions from the ion trap is controlled by electrodes 13 and 14.

Collision cell 16 can be any of a variety of means to fragment parent ions into daughter ions. Preferably, the fragmentation cell used keeps the ions contained along a path leading to the second stage analyzer, typically a TOF analyzer, downstream. As shown, the collision cell 16 typically has a RF only multi-pole 17 therein. Ions are thus focused in center region 36 and make collision with Argon or other collision gas in the cell. This process, providing one non-exclusive form of ion fragmenter is referred to as a collision induced dissociation (CID) device. The daughter ions passing out of the fragmenter through an exit (also called fragment ions or product ions) are then typically focused and cooled by another RF only multi-pole ion guide 19 and preferably pass through an electrostatic lens and ion gate assembly 20, 21, 22 before entering in input into the second stage mass analyzer, preferably in the form of a time of flight (TOF) device.

In the TOF device, a push pulse (i.e. 300V) is applied on electrode 23. Ions are pushed to the acceleration region 25. The potential difference between mesh 26 and 24 accelerates ions to high speed. Ions will fly at a constant speed through a field-free drift region 27, and then are reflected by a reflectron, also called an ion mirror 28-30, before finally striking onto a multi channel plate 32 (MCP) or other detector. Ion striking signals are typically detected by an anode 33 located behind the MCP. The pusher pulse of the acceleration region 25, typically 10Hz to 20kHz, also triggers a

timing reference for a digitizer. Based on time difference between ion arriving signal and reference trigger signal, time of ion flight is recorded digitally into a computer, later to be converted to mass/charge ratio data for that ion. The computer is configured as one form of a means to acquire, organize, store and/or display the data as depicted in FIGS. 4-10.

- 5 The TOF mass analyzer beneficially very quickly scans the daughter ions so that the TOF device is ready to scan daughter ions from the next parent ion subsequently entering the collision cell. To keep the overall tandem mass spectrometer functioning properly in real time, the TOF device preferably scans at least one hundred times faster than the first stage mass analyzer, and preferably one thousand times faster. The first stage mass analyzer can be in the form of a slow
- 10 TOF device with a gate style detector that can pass parent ions to the collision cell before resulting daughter ions are analyzed by a first TOF device, to keep the speed differential between the two mass analyzers sufficient to avoid overlap of daughter ions from different parent ions in the second stage TOF device.

- In FIG. 2, infrared multi-photon dissociation (IRMPD) is used to fragment ions. A laser beam
- 15 from an infrared laser 39 is reflected by a mirror 40 into a RF only multi-pole region. A parent ion beam from the ion trap is deflected by a deflector 41 into the same region. The parent ions are fragmented by IR radiation. Unlike CID, IRMPD does not require certain kinetic energy for parent ions, and does not need collision gas. Otherwise, the tandem mass spectrometer embodiment of FIG. 2 is similar to that of FIG. 1. The fragmenter can similarly be designed to
- 20 operate on the principals of collisionally activated dissociation or surface induced dissociation, to achieve the dividing of the parent ions into daughter ions.

- FIG. 3 is a timing diagram for the ion trap-TOF tandem MS/MS apparatus of this invention depicting time advance from left to right. Ion gate 9 (FIG. 1) drops the voltage 50 (FIG. 3) to allow ions to enter into the ion trap 37. The ion gate 9 (FIG. 1) will stay open for a short amount
- 25 of time 52 (i.e. 1-15 milliseconds), then it will close by rising the voltage 51 (FIG. 3). Meanwhile, the ion trap will trap ions 54. During the ion trap scan cycle 53 that follows, the ion gate 9 will remain closed 59. Ions upstream of the ion gate 9 can be accumulated in a multi-pole ion guide with an electric gate if desired, as described above, but are kept out of the ion trap 37. Also, during the ion trap scan cycle 53, voltage pulses 55, 56 of typically approximately 300V will be sent to the
- 30 TOF pusher 23 (FIG. 1) to start the TOF scan.

- Each TOF scan represents ions ejected from the ion trap between the last pulse and current pulse, which is a small slice 57 of parent ions. Additional pulses will result in additional slices of parent scans 58 for different parent ion mass/charge ratios. The resolution of such slices will depend on ion trap scan speed and TOF pusher pulse frequency. For example, if trap scan rate is
- 35 2000 amu/sec, that will scan from 300-1300 in half a second, and if TOF pusher frequency is 20kHz, that will give 0.1 amu resolution for parent ions. There will be a few microseconds time delay to allow ions through the collision cell, and also have some velocity variations during this transition, affecting parent ion resolution slightly. If IRMPD fragmentation is used, the daughter



ions remain closer together and parent ion resolution is not so affected.

FIG. 4 is a three-dimensional fragmentation spectrum of a five angiotensons mixture sample as it would appear if analyzed using the tandem mass spectrometer of this invention. The x-axis 60 represents daughter ion (fragmentation ion) mass to charge ratios and the y-axis 61 represents parent ion mass to charge ratios. The data shown is actually compiled from multiple separate analyses with prior art apparatuses and combined in a fashion depicting how this invention would collect and display data in a single analysis. FIG. 4 graphically illustrates the multi-dimensional information of a complete parent-daughter MS/MS map. The spectrum in FIG. 4 shows five peptides with different adducts, also charge states are from one to three. This spectrum represents the complexity with which multiple compounds may co-elute from a single HPLC peak. It only takes a few seconds by this invention to get a complete 2-D spectrum as shown. In contrast, hours of extensive scanning would be required with prior art tandem mass spectrometry.

FIG. 5 is a two-dimensional plot from the same spectrum of data shown in FIG. 4. Once data from this spectrum has been entered into a digital computer it can be viewed in various ways to provide the desired information. For instance, if the data along a horizontal line 66 is plotted alone, as shown in FIG. 6, a daughter scan spectrum for parent ion  $m/z=884$  is provided. If the data along a vertical line 65 (FIG. 5) is plotted alone, as shown in FIG. 7, a parent scan for daughter ion  $m/z=610$  is provided.

A diagonal line 67 with the x coordinate equal to the y coordinate represents the data related to unfragmented parent ions. A diagonal line 68 to the left of this first diagonal line 67 where the x coordinate is 18 less than the y coordinate, represents what a neutral loss scan plotted alone would provide, as shown in FIG. 8. If the data for the sum of all x values is plotted on the y-axis, as shown in FIG. 10, a full MS scan of the parent ions is provided.

FIG. 9 shows the 2-D neutral loss map from the same spectrum. In this plot, the data points are shifted to the left. The distance of the shifting is equal to the y value. As a result, the new x-axis 85 is transformed to the value of neutral loss. Line 68 in FIG. 5 becomes line 88 in this plot of FIG. 9. Line 67 (FIG. 5) becomes line 89 (FIG. 9). This plot of FIG. 9 gives a clear two-dimensional picture that graphically illustrates the neutral loss relations of each parent ion. Every point lined up vertically (i.e. at 90) represents the same neutral loss.

This disclosure is provided to reveal a preferred embodiment of the invention and a best mode for practicing the invention. Having thus described the invention in this way, it should be apparent that various different modifications can be made to the preferred embodiment without departing from the scope and spirit of this disclosure. When structures are identified as a means to perform a function, the identification is intended to include all structures which can perform the function specified. When structures of this invention are identified as being coupled together, such language should be interpreted broadly to include the structures being coupled directly together or coupled together through intervening structures. Such coupling could be permanent or temporary and either in a rigid fashion or in a fashion which allows pivoting, sliding or other relative motion

while still providing some form of attachment. When elements are described as upstream or downstream relative to other elements, the elements can be directly upstream or downstream with no intervening elements or indirectly upstream or downstream with intervening elements therebetween.

5

### **Industrial Applicability**

This invention exhibits industrial applicability in that it provides apparatuses and methods for more rapidly, more completely, more flexibly and more efficiently obtaining data of the type obtained by tandem mass spectrometry (MS/MS).

Another object of the present invention is to provide an apparatus and method for rapidly obtaining ion mass data with high sensitivity and a large dynamic range.

Another object of the present invention is to provide a single mass spectrometry instrument that has good versatility and can perform in multiple scan modes.

Another object of the present invention is to provide a method and apparatus for obtaining MS/MS type two dimensional data about parent ions and daughter ions sufficiently rapidly to facilitate combination with a chromatographic apparatus, such that complete multidimensional data can be obtained in real time, during the relatively short duration of a single chromatographic peak.

Another object of the present invention is to provide a method and apparatus that uses a non-destructive mass analyzer for both first and second stage analysis for obtaining complete spectrum MS/MS type data.

Other further objects of this invention, which demonstrate its industrial applicability, will become apparent from a careful reading of the included detailed description, from a review of the enclosed drawings and from review of the claims included herein.

30

35

**CLAIMS**

What is claimed is:

- 5      Claim 1 - A multiple stage mass spectrometer, comprising in combination:  
         a first non-destructive mass analyzer including an ion trap;  
         an ion fragmenter downstream from the first mass analyzer; and  
         a second mass analyzer downstream from said ion fragmenter.
- 10      Claim 2 - The mass spectrometer of Claim 1 wherein said ion trap includes an inlet coupled to  
a source of a sample including at least one species of parent ions.
- Claim 3 - The mass spectrometer of Claim 2 wherein said ion trap is a three dimensional ion  
trap.
- 15      Claim 4 - The mass spectrometer of Claim 2 wherein said ion trap is a linear ion trap.
- Claim 5 - The mass spectrometer of Claim 1 wherein said first mass analyzer includes an  
outlet for the parent ions contained therein, said outlet aligned with said ion fragmenter.
- 20      Claim 6 - The mass spectrometer of Claim 5 wherein said ion fragmenter is adapted to divide  
parent ions from said first mass analyzer into daughter ions.
- Claim 7 - The mass spectrometer of Claim 6 wherein said ion fragmenter includes a collision  
25 induced dissociation device.
- Claim 8 - The mass spectrometer of Claim 6 wherein said fragmenter includes an infrared  
multi-photon dissociation device.
- 30      Claim 9 - The mass spectrometer of Claim 6 wherein said fragmenter includes a collisionally  
activated dissociation device.
- Claim 10 - The mass spectrometer of Claim 6 wherein said ion trap includes an outlet aligned  
with an entrance into said ion fragmenter.
- 35      Claim 11 - The mass spectrometer of Claim 6 wherein said fragmenter includes an exit aligned  
with said second mass analyzer.

Claim 12 - The mass spectrometer of Claim 6 wherein said second mass analyzer is adapted to separate ions on the order of at least one thousand times faster than separation by said first mass analyzer.

5      Claim 13 - The mass spectrometer of Claim 12 wherein said second mass analyzer includes a time of flight device.

10      Claim 14 - The mass spectrometer of Claim 1 wherein a computation device, a memory and a display are coupled to both said first mass analyzer and said second mass analyzer to receive data from said first mass analyzer and said second mass analyzer associated with ions detected by said first mass analyzer and said second mass analyzer.

15      Claim 15 - The mass spectrometer of Claim 14 wherein said computation device is adapted to combine data from said first mass analyzer with data from said second mass analyzer to create a two dimensional output plot of mass to charge ratios for at least two parent ions detected by said first mass analyzer and at least two daughter ions detected by said second mass analyzer.

20      Claim 16 - The mass spectrometer of Claim 1 wherein a source of ions is provided upstream from an inlet into said first mass analyzer, said source of ions including an electrospray upstream from a RF only multi-pole ion guide upstream of an electrostatic lens upstream of an inlet into said first mass analyzer.

25      Claim 17 - A tandem mass spectrometer, comprising in combination:  
a first non-destructive mass analyzer having a parent ion inlet and a parent ion outlet;  
means to fragment parent ions downstream from said first mass analyzer into daughter ions; and  
a second mass analyzer having a daughter ion input downstream from said parent ion fragmenting means.

30      Claim 18 - The mass spectrometer of Claim 17 wherein said first mass analyzer outlet is selectively openable and closable with said first mass analyzer adapted to release parent ions therefrom when said ion outlet is open and to retain parent ions within said first mass analyzer when said outlet is closed.

35      Claim 19 - The mass spectrometer of Claim 18 wherein said first mass analyzer includes an ion trap.

Claim 20 - The mass spectrometer of Claim 19 wherein said ion trap is a three dimensional ion trap.

Claim 21 - The mass spectrometer of Claim 19 wherein said ion trap is a linear ion trap.

5

Claim 22 - The mass spectrometer of Claim 19 wherein said parent ion outlet of said first mass analyzer is adjusted from an open position to a closed position by adjusting a voltage of an electric field of said ion trap.

10

Claim 23 - The mass spectrometer of Claim 17 wherein both said first mass analyzer and said second mass analyzer are coupled to a means to acquire data related to mass to charge ratios of ions detected by both of said mass analyzers, with said data displayed in two dimensions.

15

Claim 24 - The mass spectrometer of Claim 23 wherein said means to acquire data includes mean to display said data in two dimensions including an x axis and a y axis with one of said axes representing a mass to charge ratio of parent ions and the other of said axes representing a mass to charge ratio of daughter ions.

20

Claim 25 - The mass spectrometer of Claim 17 wherein said second mass analyzer includes a time of flight device.

Claim 26 - The mass spectrometer of Claim 25 wherein said first mass analyzer includes an ion trap.

25

Claim 27 - The mass spectrometer of Claim 17 wherein said fragmentation means includes a collision cell.

30

Claim 28 - The mass spectrometer of Claim 26 wherein said collision cell includes an RF only multi-pole therein and a collision gas therein such that said collision cell is adapted to cause collision induced dissociation.

35

Claim 29 - The mass spectrometer of Claim 17 wherein said fragmentation means includes an infrared laser oriented to expose the parent ions with photons to induce fragmentation of the parent ions into daughter ions.

Claim 30 - The mass spectrometer of Claim 16 wherein said first mass analyzer includes a time of flight device, said second mass analyzer includes a time of flight device separate from said first mass analyzer, said second mass analyzer configured to operate at least one hundred times faster

than said second mass analyzer.

Claim 31 - A two stage mass analyzer, comprising in combination:

5 a first mass analyzer in the form of an ion trap, said ion trap having an inlet downstream from a parent ion source and a parent ion outlet;

a fragmentation cell downstream from said ion trap, said fragmentation cell adapted to divide parent ions from said ion trap outlet into daughter ions, said fragmentation cell including a daughter ion exit; and

10 a second mass analyzer including a time of flight device downstream from said fragmentation cell exit.

Claim 32 - The two stage mass analyzer of Claim 31 wherein a computer is coupled to said first mass analyzer and said second mass analyzer, said computer adapted to acquire mass to charge ratio data for both parent ions from said first mass analyzer and daughter ions from said  
15 second mass analyzer.

Claim 33 - The two stage mass analyzer of Claim 32 wherein said computer is adapted to correlate parent ion data with daughter ion data.

20 Claim 34 - The two stage mass analyzer of Claim 33 wherein said computer is adapted to display correlated parent ion and daughter ion data in the form of a two dimensional plot including an x axis and a y axis with one of said axes representing a mass to charge ratio of the parent ions and the other of said axes representing the mass to charge ratio of the daughter ions.

25 Claim 35 - The two stage mass analyzer of Claim 31 wherein said fragmentation cell includes a collision induced dissociation device.

Claim 36 - The two stage mass analyzer of Claim 31 wherein said fragmentation cell includes an infrared multi-photon dissociation device.  
30

Claim 37 - The two stage mass analyzer of Claim 31 wherein said fragmentation cell includes a collisionally activated dissociation device.

Claim 38 - The two stage mass analyzer of Claim 31 wherein said ion source upstream of said  
35 first mass analyzer is an output of a chromatography device.

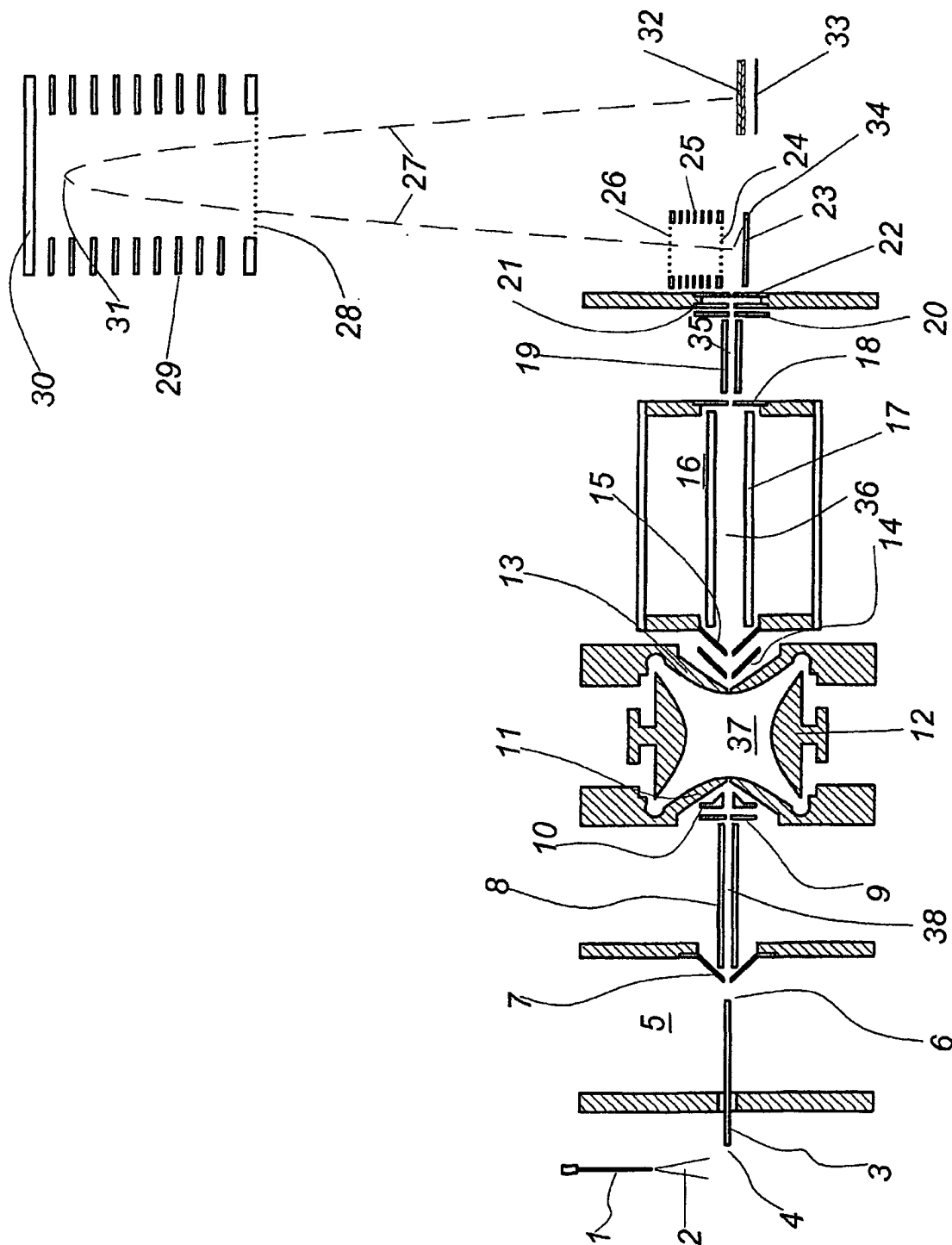


Figure 1

2/10

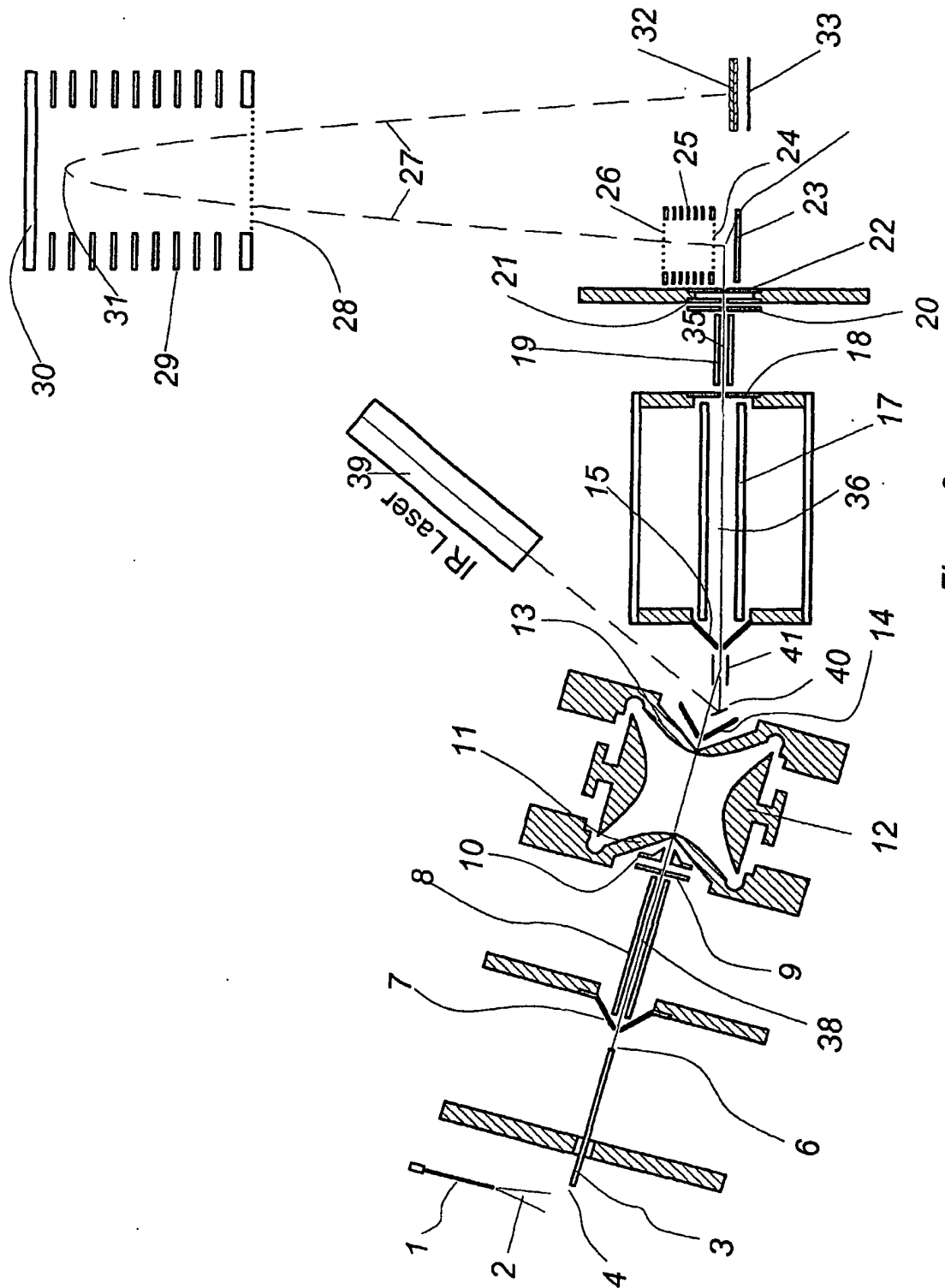


Figure 2



3/10

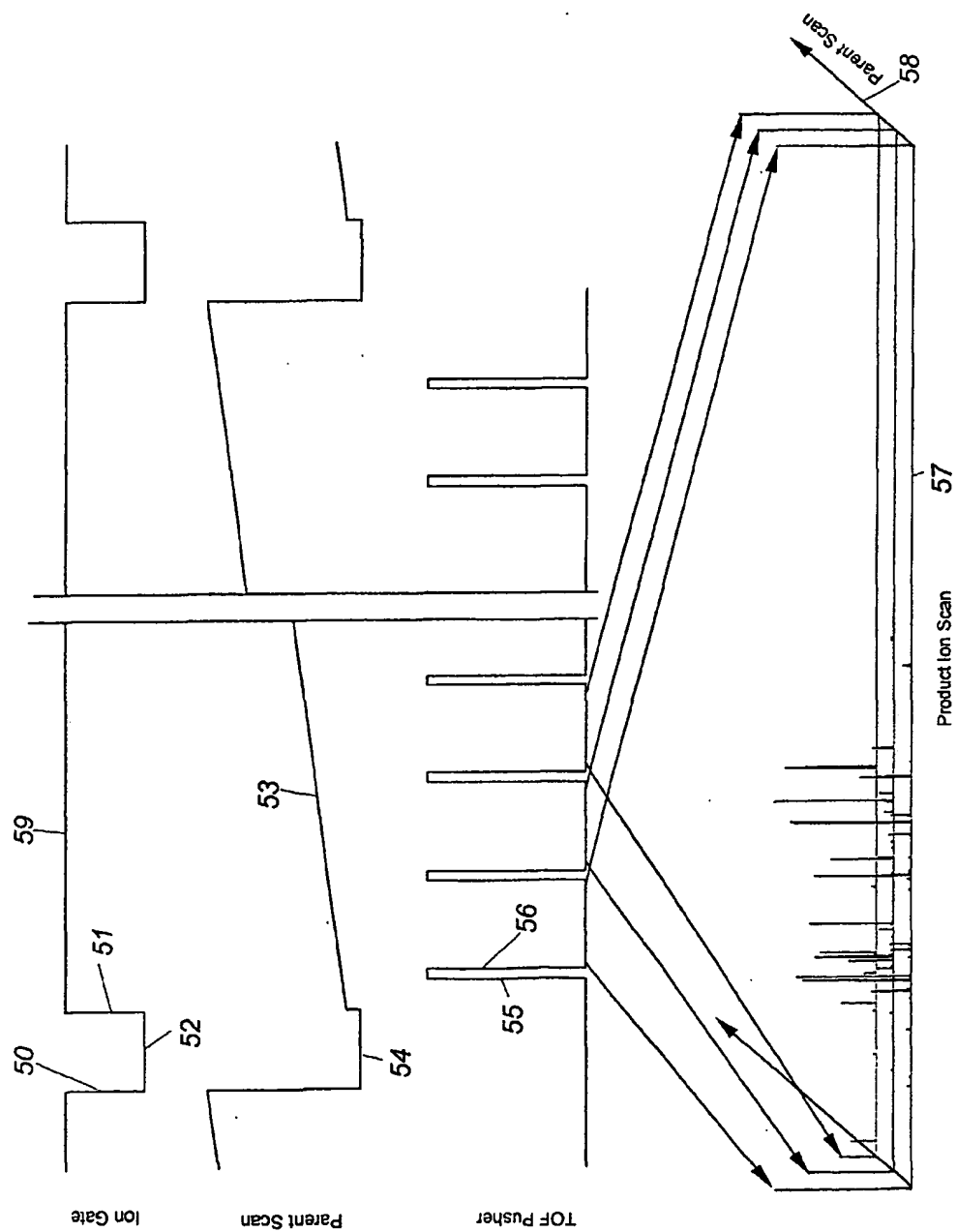


Figure 3

4/10

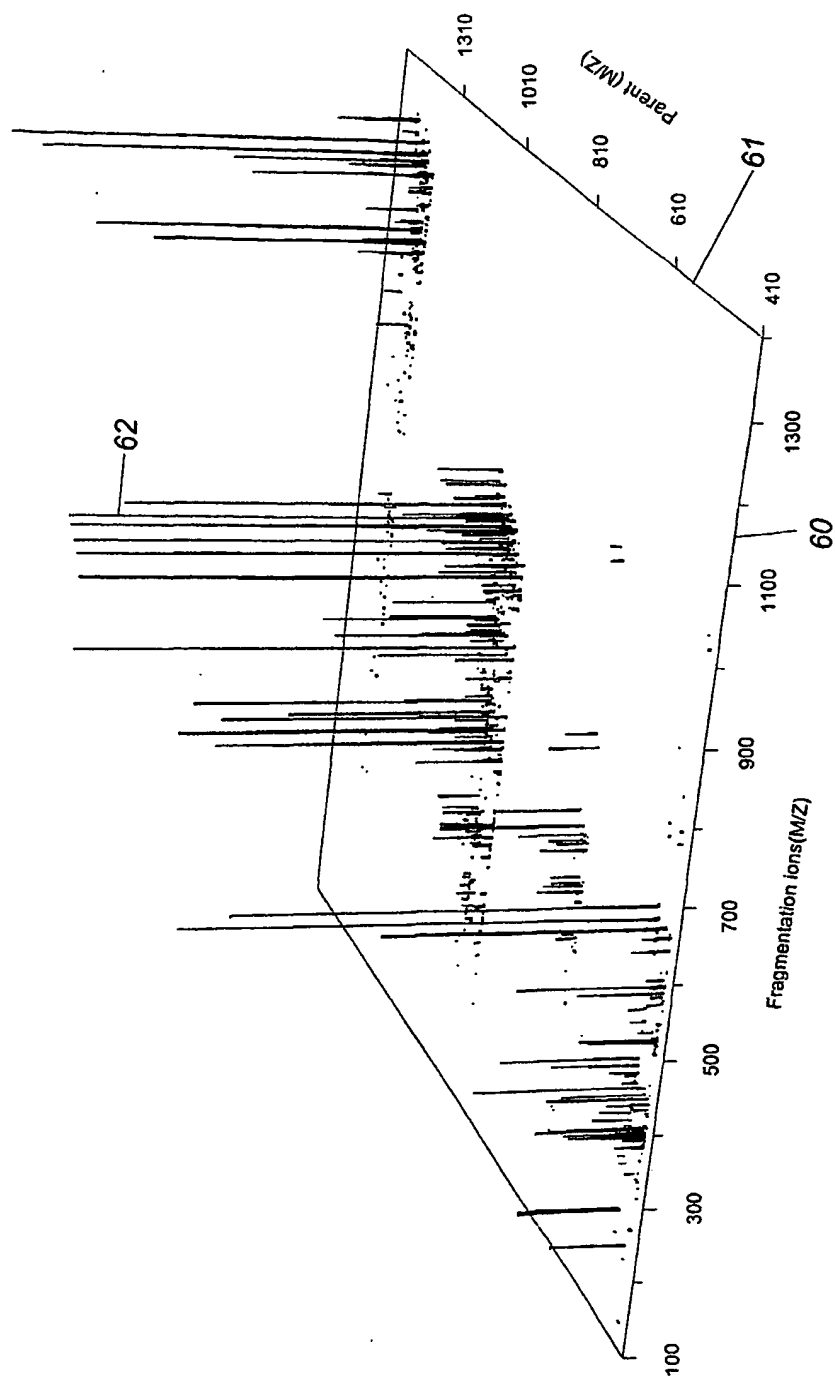


Figure 4

5/10


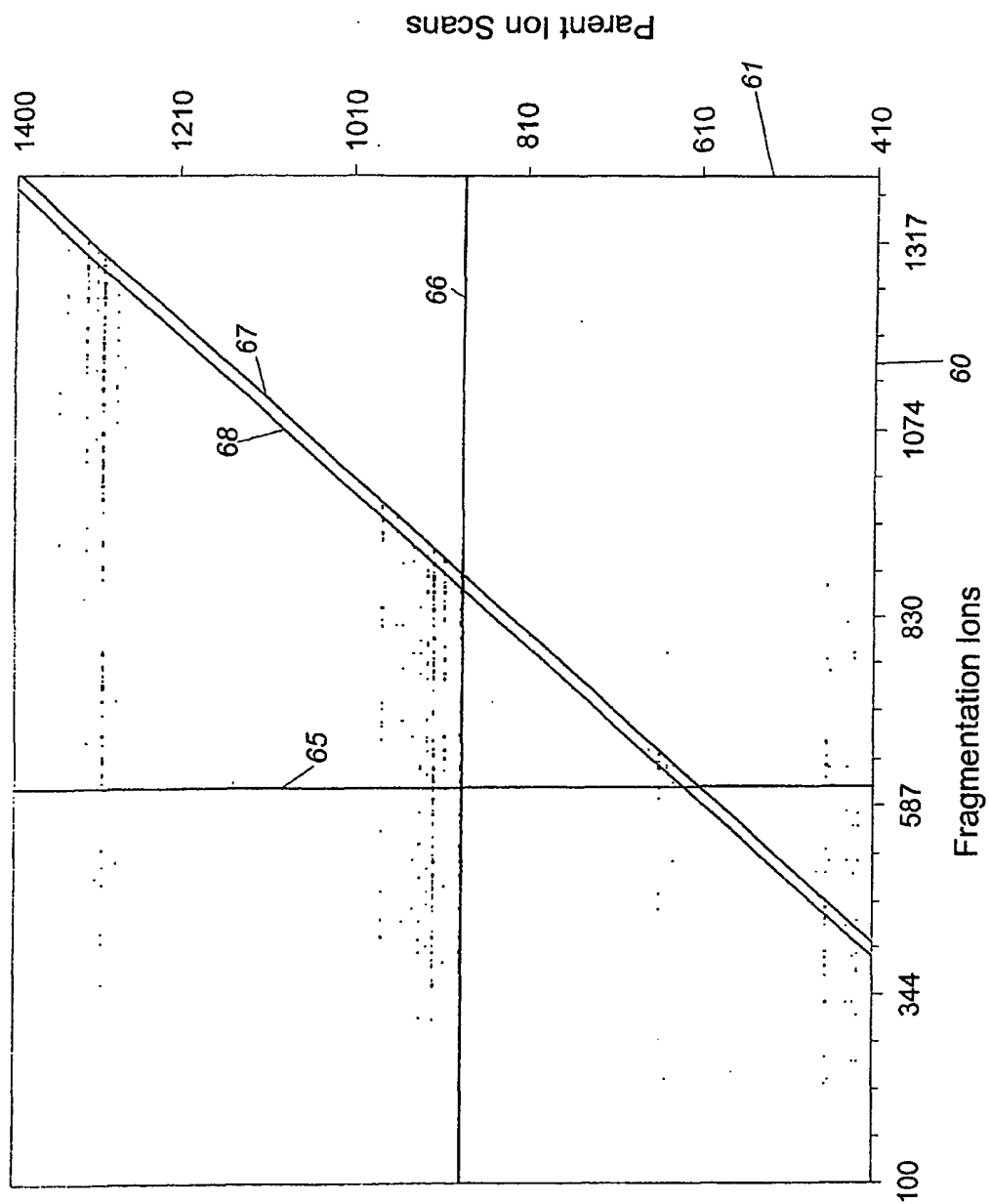


Figure 5

6/10

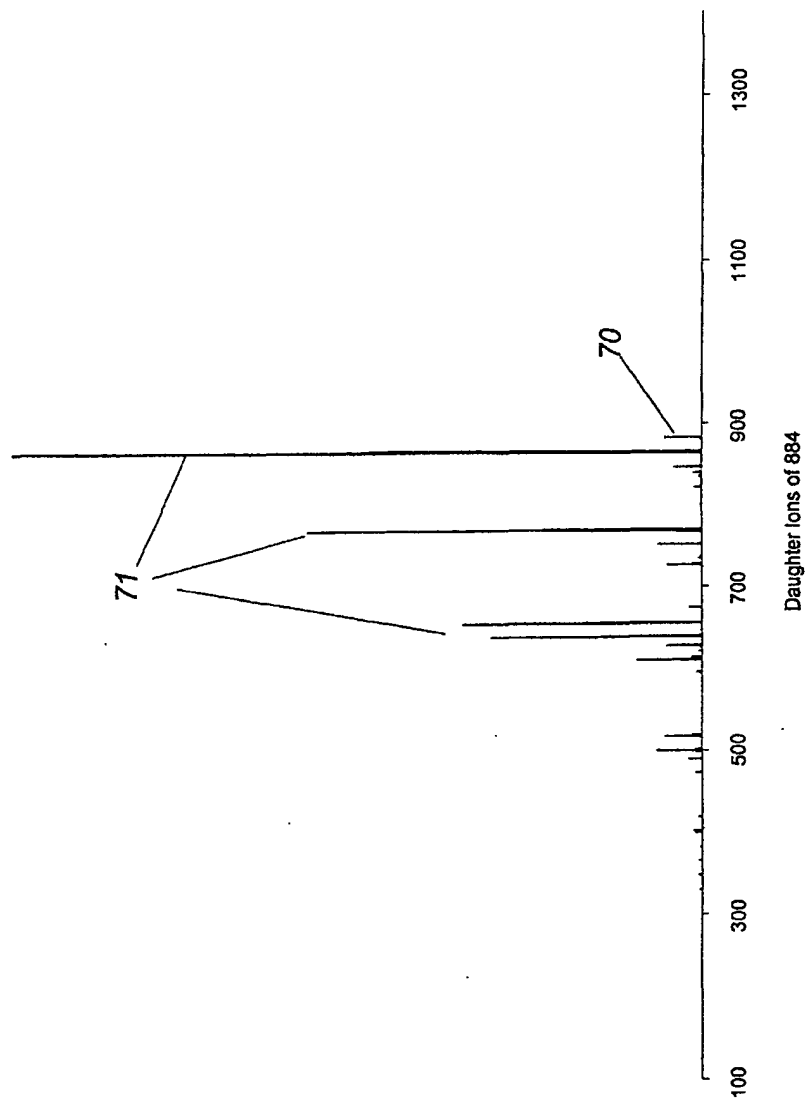

7/10

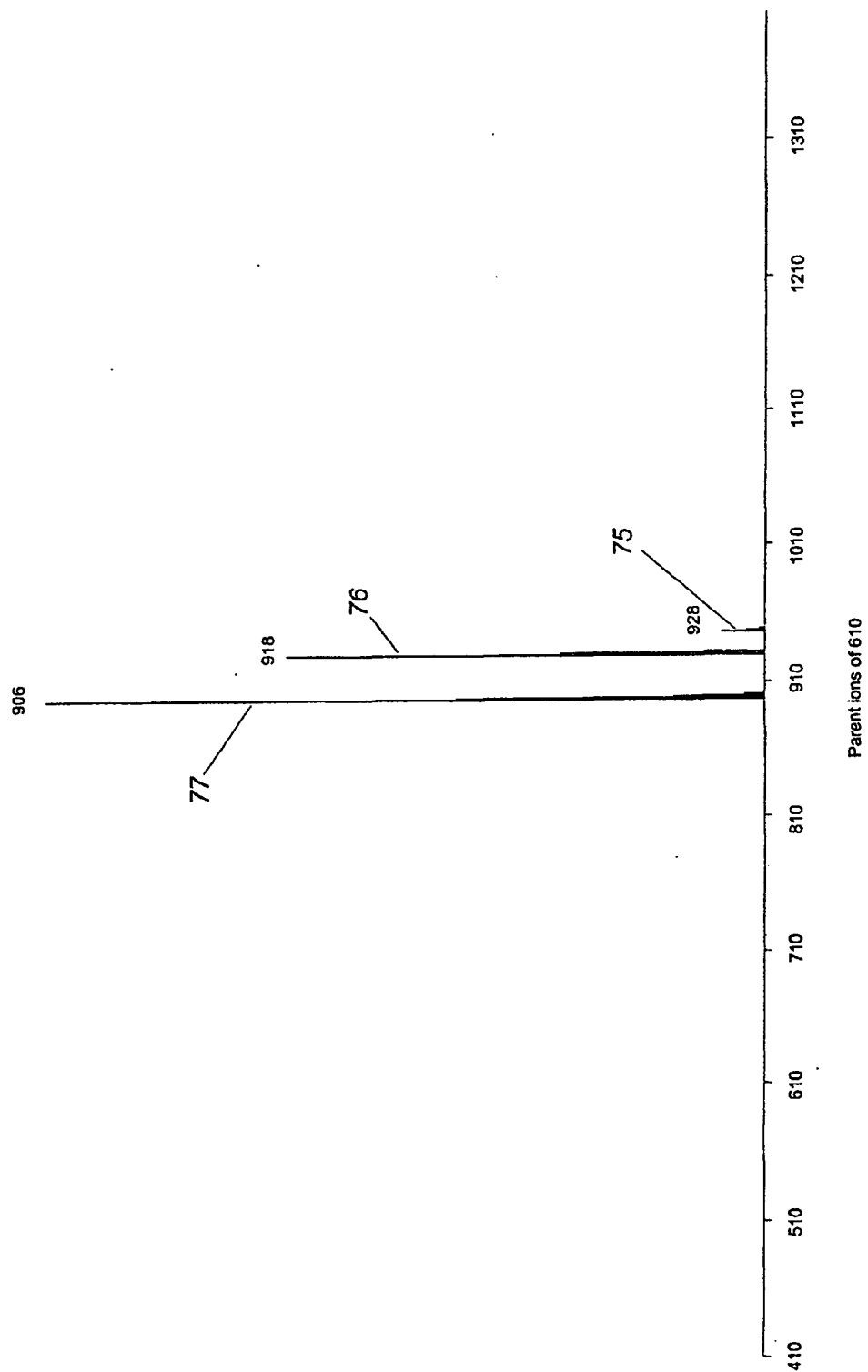


Figure 7

8/10

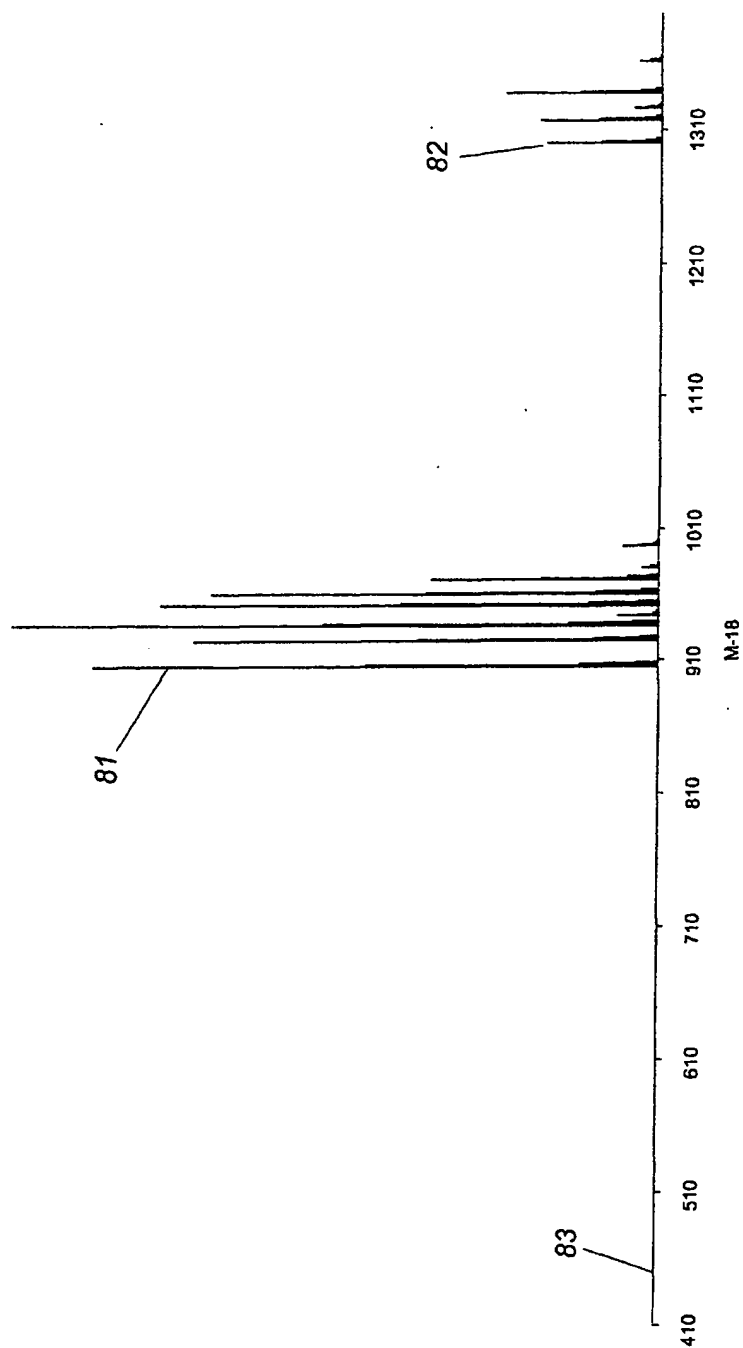


Figure 8

9/10

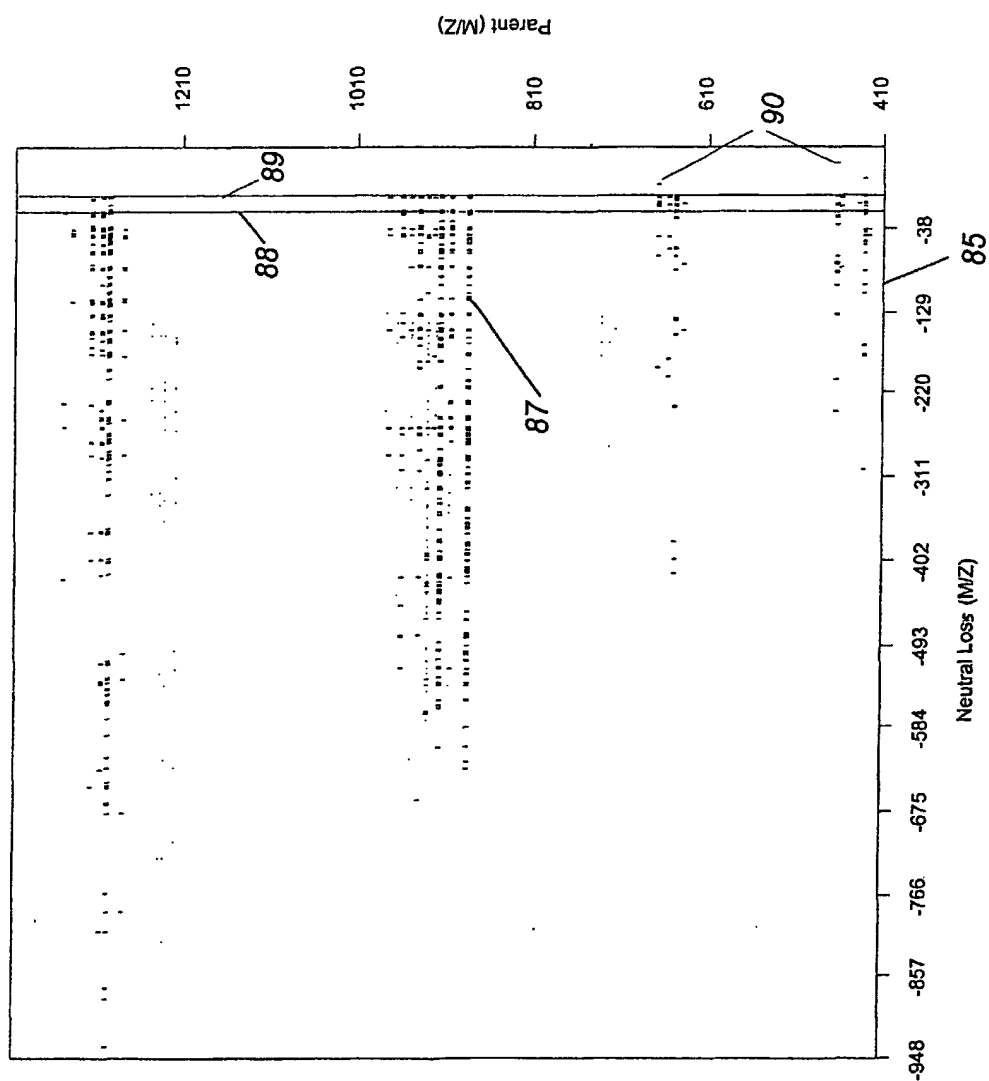


Figure 9

10/10

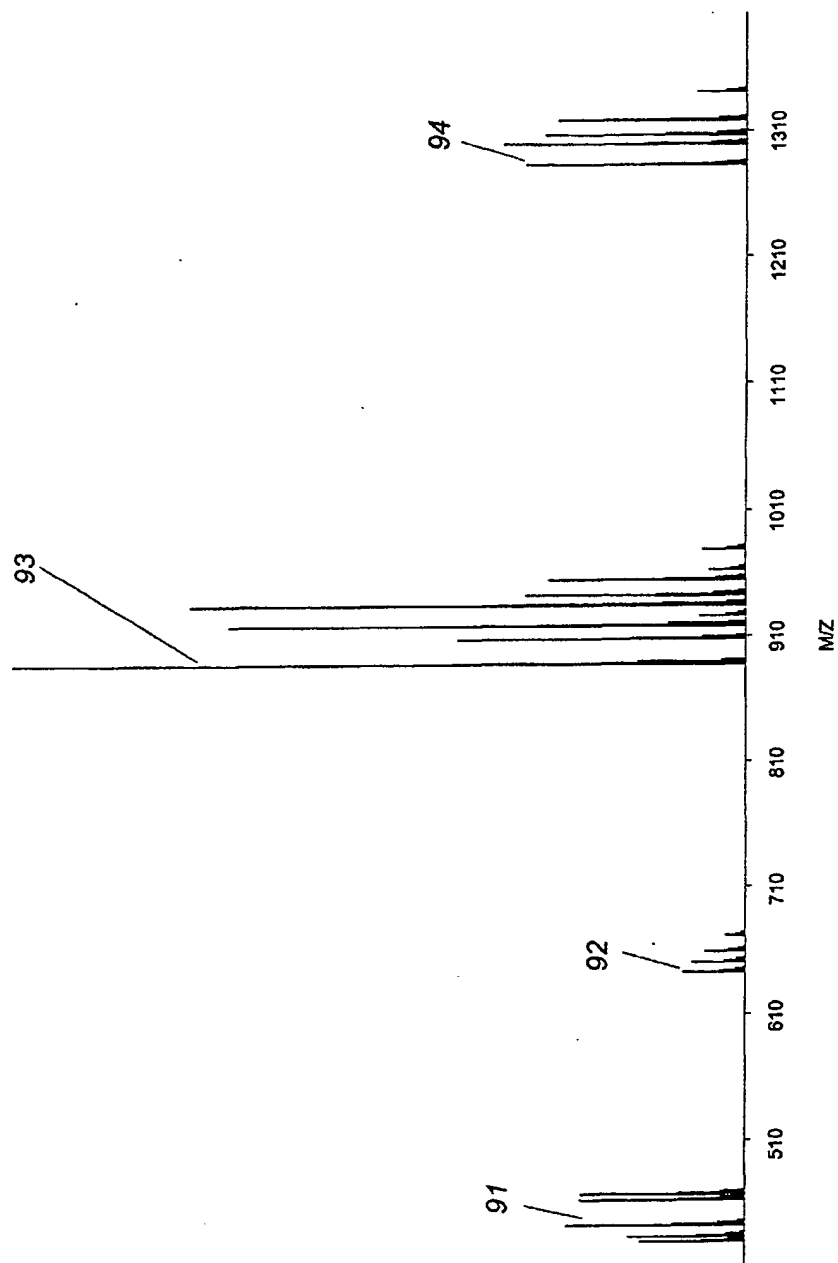


Figure 10



## INTERNATIONAL SEARCH REPORT

PCT/US 03/15718

**A. CLASSIFICATION OF SUBJECT MATTER**  
IPC 7 H01J49/42 H01J49/40

According to International Patent Classification (IPC) or to both national classification and IPC

**B. FIELDS SEARCHED**

Minimum documentation searched (classification system followed by classification symbols)

IPC 7 H01J

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the International search (name of data base and, where practical, search terms used)

EPO-Internal, WPI Data, PAJ, INSPEC

**C. DOCUMENTS CONSIDERED TO BE RELEVANT**

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y	EP 0 898 297 A (MICROMASS LTD) 24 February 1999 (1999-02-24)	1-7, 9-11, 17-22, 24,25, 27,28, 31,35, 37,38
A	column 9 -column 10  --- -/--	8,12-15, 23,29,30



Further documents are listed in the continuation of box C.



Patent family members are listed in annex.

## \* Special categories of cited documents:

- \*A\* document defining the general state of the art which is not considered to be of particular relevance
- \*E\* earlier document but published on or after the international filing date
- \*L\* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- \*O\* document referring to an oral disclosure, use, exhibition or other means
- \*P\* document published prior to the international filing date but later than the priority date claimed

- \*T\* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- \*X\* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- \*Y\* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- \*G\* document member of the same patent family

Date of the actual completion of the international search

2 September 2003

Date of mailing of the international search report

09/09/2003

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2  
NL - 2280 HV Rijswijk  
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl.  
Fax: (+31-70) 340-3016

Authorized officer

Hulne, S

## INTERNATIONAL SEARCH REPORT

PCT/US 03/15718

## C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y	WILHELM U ET AL: "Ion storage combined with reflectron time-of-flight mass spectrometry: ion cloud motions as a result of jet-cooled molecules" INTERNATIONAL JOURNAL OF MASS SPECTROMETRY AND ION PROCESSES, ELSEVIER SCIENTIFIC PUBLISHING CO. AMSTERDAM, NL, vol. 152, no. 2, 29 February 1996 (1996-02-29), pages 111-120, XP004036547 ISSN: 0168-1176 page 112 -page 113 ----	1-7, 9-11, 17-22, 24,25, 27,28, 31,35, 37,38
A	GOERINGER E: "Tandem quadrupole time-of-flight instrument for mass spectrometry / mass spectrometry" ANALYTICAL CHEMISTRY, AMERICAN CHEMICAL SOCIETY. COLUMBUS, US, vol. 56, 1984, pages 2291-2295, XP002151968 ISSN: 0003-2700 figures 1,2 ----	1,17,31
A	MICHAEL S M ET AL: "AN ION TRAP STORAGE/TIME-OF-FLIGHT MASS SPECTROMETER" REVIEW OF SCIENTIFIC INSTRUMENTS, AMERICAN INSTITUTE OF PHYSICS. NEW YORK, US, vol. 63, no. 10 PT 1, 1 October 1992 (1992-10-01), pages 4277-4284, XP000316284 ISSN: 0034-6748 page 4278 -page 4279; figures 3,4 -----	1,17,34

## INTERNATIONAL SEARCH REPORT

PCT/US 03/15718

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
EP 0898297 A	24-02-1999	CA 2245465 A1	22-02-1999
		EP 0898297 A2	24-02-1999
		JP 11154486 A	08-06-1999
		US 6107623 A	22-08-2000
<hr/>			